



FAQ

AED PAD PLACEMENT

Frequently asked questions when using AED Electrode Pads.

Where do you place the electrode pads?


This depends on the age and size of the victim.

For adults place the first pad on the upper right-hand side of the chest below the collar bone.

Place the second pad on the lower left side of the chest (below the pectoral muscle of a male patient or under the breast of a female patient)

ADULTS

(includes children
8+ years or more
than 25 kg)

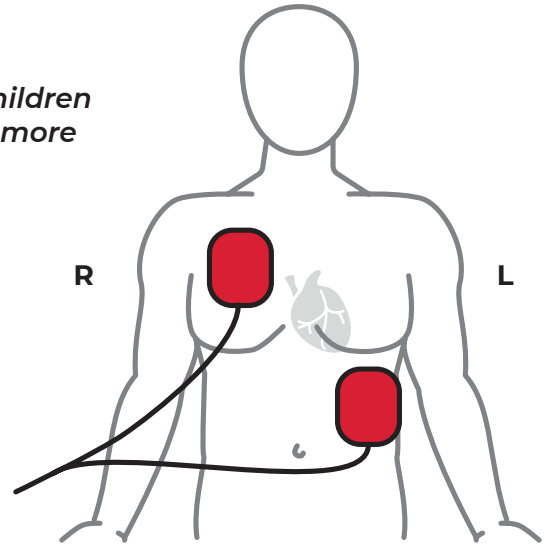

Can I use an AED on a child / baby?

Yes, some AEDs have paediatric pads which are smaller and can be placed as per instructions. If there are no paediatric pads available you can use the adult pads on larger children. For smaller children and babies you can place these pads front and back, keeping the heart in a direct line between the pads.

CHILDREN


(less than 8 years
or less than 25 kg)

USING ADULT PADS

For smaller children and babies place these pads front and back, keeping the heart in a direct line between the pads.


USING PAEDIATRIC PADS

Place as per instructions on the pads.

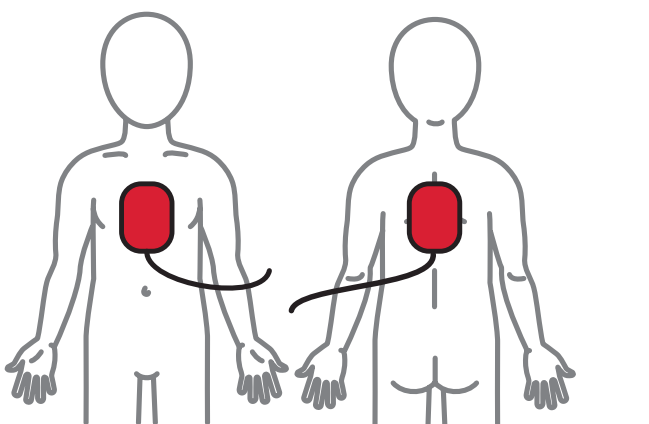

FRONT

BACK

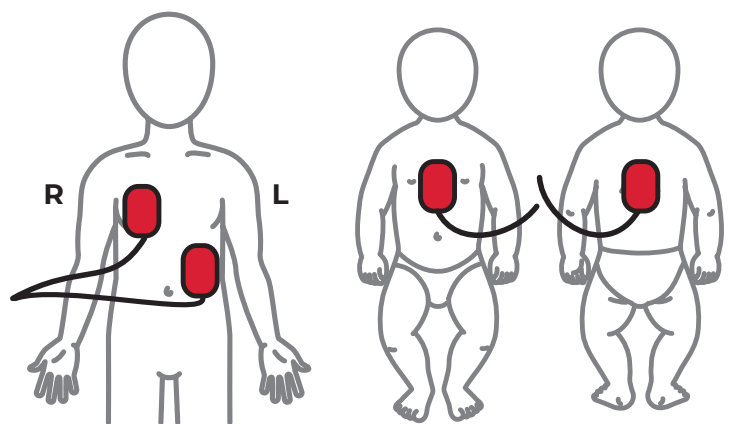

FRONT

BACK

0800 233 342

www.heartsaver.co.nz